---

# CERTEST

## *Crypto+Giardia* *+Entamoeba*

---

ONE STEP

*Cryptosporidium, Giardia and*  
*Entamoeba*

COMBO CARD TEST

---

CERTEST BIOTEC S.L.

---

## CERTEST *Crypto+Giardia+Entamoeba*

One Step test to detect *Cryptosporidium*, *Giardia* and *Entamoeba* in combo card format /  
Prueba combo de un solo paso para detección de *Cryptosporidium*, *Giardia* y *Entamoeba* en formato cassette

### ENGLISH

#### INTENDED USE

CerTest *Crypto+Giardia+Entamoeba* one step combo card test is a coloured chromatographic immunoassay for the simultaneous qualitative detection of *Cryptosporidium*, *Giardia* ( $\alpha$ 1-*giardin* and/or CWP1), *Entamoeba* (*histolytica* and/or *dispar*) in stool samples.

CerTest *Crypto+Giardia+Entamoeba* combo card test offers a simple and highly sensitive screening assay to make a presumptive diagnosis of cryptosporidiosis and/or giardiasis and/or amoebiasis.

#### INTRODUCTION

*Cryptosporidium parvum*, *Giardia lamblia* and *Entamoeba histolytica* are three of the major causes of protozoan-induced diarrheal disease and are the most frequently identified protozoan parasites causing waterborne disease outbreaks. *Giardia* causes an intestinal illness called giardiasis. *Cryptosporidium* is responsible for a similar illness called cryptosporidiosis. And amoebiasis is the infection of the human gastrointestinal tract by *Entamoeba histolytica*. These infections have become the most common causes of waterborne diseases (found in both drinking and recreational water) in humans.

*Giardia*, a flagellated protozoan, inhabits the upper part of the small intestine of its host and has a two major states in the life cycle: trophozoites which produces the antigens ( $\alpha$ -1 *giardin*) and cyst with produces the antigens (CWP1) After the host ingest the cysts, which are the infective stage, the trophozoites emerge from the cysts in the duodenum and attach to the small intestinal mucosa. They undergo mitotic division in the intracellular lumen, some will encyst to protect themselves and will be eliminated from the host in the feces. The trophozoite is the vegetative form and replicates in the small intestine.

Giardiasis is a diarrheal illness caused by a very small parasite, *Giardia intestinalis* (also known as *Giardia lamblia* and *Giardia duodenalis*). Once an animal or person is infected with *Giardia*, the parasite lives in the intestine and is passed in the stool. The parasite is protected by an outer shell and can survive outside the body and in the environment for a long time. The most common symptoms of giardiasis include: diarrhea, loose or watery stool, stomach cramps and upset stomach. These symptoms generally begin 1-2 weeks after infection and may last 2-6 weeks in healthy individuals. Sometimes symptoms last longer and may lead to weight loss and dehydration. Some people will have no symptoms. However, people with weakened immune systems (e.g., persons with HIV/AIDS, cancer patients and transplant patients) or the elderly may have a more serious infection that can lead to severe illness or death.

Several species and genotypes of *Cryptosporidium* have been associated with human infections, although *C. hominis* (previously known as the human genotype of *C. parvum*) and *C. parvum* have been identified as the chief causes of human cryptosporidiosis. The disease is more common in young children and immunosuppressed patients, in whom it can be very serious, but it can also affect healthy people of any age. However, in immunocompetent patient, clinical manifestations of cryptosporidiosis range from the absence of symptoms, though mild symptoms, to a severe gastroenteritis-like syndrome with watery diarrhoea. *Cryptosporidium parvum* is the major cause of persistent diarrhoea in developing countries. This parasite is recognised as a highly infectious enteric pathogen and infective stage is transmitted by the fecal-oral route. Symptoms of cryptosporidiosis include watery diarrhoea, stomach cramps, weight loss, nausea and sometimes fever. The incubation period is 5 to 28 days with a mean of 7.2 days. Other symptoms include abdominal pain and fatigue.

*Entamoeba histolytica*, *E. dispar*, and *E. moshkovskii* are morphologically identical but biochemically and genetically different. These parasites colonize the human gut, but only *E. histolytica* is thought to be capable of causing disease.

Amoebiasis is the infection of the human gastrointestinal tract by *Entamoeba histolytica*, a protozoan parasite that is capable of invading the intestinal mucosa and may spread to other organs, mainly the liver. The acceptance of *E. dispar* as a distinct but closely related protozoan species has had profound implications for the epidemiology of amoebiasis, since most asymptomatic infections found worldwide are now attributed to this non invasive amoeba. Roughly, 90% of people infected with *Entamoeba* infections have *E. dispar*.

The disease may manifest itself as an acute, chronic or as an asymptomatic infection. The leading symptom is dysentery whose main consequence one galling of the colon, diarrhoea and abdominal pain. Complications of the disease may be hepatic abscesses, pulmonary abscesses or even cerebral abscesses which, if untreated, usually end in death.

Infection occurs by water pollution, vegetables, raw fruits and other food or poorly cooked or washed with infective cysts from contaminated faeces. It is possible that flies and cockroaches carrying cysts from faeces to food

Amoebic dysentery occurs frequently in tropical countries but are also cases in temperate and cold areas. In Africa, tropical Asia and Latin America, more than two thirds of the population has these intestinal parasites, although most can be virtually asymptomatic infections.

#### TEST PRINCIPLE

CerTest *Crypto+Giardia+Entamoeba* is based on the principle of a qualitative immunochromatographic assay for the determination of *Cryptosporidium*, *Giardia* ( $\alpha$ 1-*giardin* and/or CWP1), *Entamoeba* (*histolytica* and/or *dispar*) in stool samples. The strip B (*Giardia*) detect trophozoites and cyst.

**Strip A** consists of a nitrocellulose membrane pre-coated with mouse monoclonal antibodies on the test line (T), in the results window, against *Cryptosporidium* and with rabbit polyclonal antibodies, on the control line (C), against a specific protein. The label/sample absorbent pad is sprayed with test label solution (mouse monoclonal antibodies anti-*Cryptosporidium*) conjugated to red polystyrene latex and control label solution (specific binding protein) conjugated to green polystyrene latex, forming two coloured conjugate complexes.

**Strip B** consists of a nitrocellulose membrane pre-coated with mouse monoclonal antibodies on the test line (T), in the results window, against *Giardia* and with rabbit polyclonal antibodies, on the control line (C), against a specific protein. The label/sample absorbent pad is sprayed with test label solution (mouse monoclonal antibodies anti-*Giardia*) conjugated to red polystyrene latex and control label solution (specific binding protein) conjugated to green polystyrene latex, forming two coloured conjugate complexes.

**Strip C** consists of a nitrocellulose membrane pre-coated with mouse monoclonal antibodies on the test line (T), in the results window, against *Entamoeba* and with rabbit polyclonal antibodies, on the control line (C), against a specific protein. The label/sample absorbent pad is sprayed with test label solution (mouse monoclonal antibodies anti-*Entamoeba*) conjugated to red polystyrene latex and control label solution (specific binding protein) conjugated to green polystyrene latex, forming two coloured conjugate complexes.

If the sample is *Cryptosporidium* positive, the antigens of the diluted sample react with the red-coloured conjugate complex (anti-*Cryptosporidium* monoclonal antibodies-red polystyrene microspheres) in the strip A, if the sample is *Giardia* positive, the antigens of the diluted sample react with the red-coloured conjugate complex (anti-*Giardia* monoclonal antibodies-red polystyrene microspheres) in the strip B and if the sample is *Entamoeba* positive, the antigens of the diluted sample react with the red-coloured conjugate complex (anti-*Entamoeba* monoclonal antibodies-red polystyrene microspheres) in the strip C, which were previously pre-dried on the absorbent pad. The mixture then moves upward on the membrane by capillary action. As the sample flows through the test membrane, the binding conjugate complexes migrate. The anti-*Cryptosporidium* antibodies present on the membrane of strip A (test line), the anti-*Giardia* antibodies present on the membrane of strip B (test line) and the anti-*Entamoeba* antibodies present on the membrane of strip C (test line) capture the coloured conjugate and the red lines will be visible in the strips. These bands are used to interpret the result.

If the sample is negative, there is no *Cryptosporidium*, *Giardia* and *Entamoeba* presence and yet, the antigens may be present in a concentration lower than the detection limit value, for which the reaction will not take place with any red-coloured conjugate complex. The anti-*Cryptosporidium*, anti-*Giardia* and anti-*Entamoeba* antibodies present on the membranes (test lines) will not capture the antigen-red-coloured conjugate complex (not formed), for which the red lines will not appear.

Whether the sample is positive or not, in the three strips, the mixture continues to move across the membranes to the immobilized specific antibodies placed in the control lines. The anti-specific protein antibodies present on three membranes will capture control green-conjugate complex and three control lines will always appear. The presence of these green lines serve as: 1) verification that sufficient volume is added, 2) that proper flow is obtained and 3) an internal control for the reagents.

### STORAGE AND STABILITY

Store as packaged in the sealed pouch at 2-30°C. The test is stable until the expiration date marked on its sealed pouch. The test must remain in the sealed pouch until use. Do not freeze.

### PRECAUTIONS

- For professional *in vitro* diagnostic use.
- Do not use after expiration date.
- All the specimens should be considered potentially hazardous and handled in the same manner as an infectious agent. A new test must be used for each sample to avoid contamination errors.
- The tests should be discarded in a proper biohazard container after testing.
- Reagents contain preservatives. Avoid any contact with skin or mucous membrane. Consult safety data sheet, available on request.
- Components provide in the kit are approved for use in the CerTest *Crypto+Giardia+Entamoeba* combo card test. Do not use any other commercial kit component.
- The presence of yellow lines in the results window (control line zone and test line zone), before using the test, is completely normal and does not imply failure of the test functionality.
- Follow Good Laboratory Practices, wear protective clothing, use disposal gloves, goggles and mask. Do not eat, drink or smoke in the working area.

### SPECIMEN COLLECTION AND PREPARATION


Stool samples should be collected in clean containers. The samples can be stored in the refrigerator (2-8°C) for 1-3 days prior to testing. For longer storage, maximum 1 month, the specimen must be kept frozen at -20°C. In this case, the sample will be totally thawed and brought to room temperature before testing. Homogenise stool sample as thoroughly as possible prior to preparation.

## CERTEST *Crypto+Giardia+Entamoeba*

*One Step test to detect Cryptosporidium, Giardia and Entamoeba in combo card format / Prueba combo de un solo paso para detección de Cryptosporidium, Giardia y Entamoeba en formato cassette*

### Specimen preparation (see illustration):

1. Take out the cap of the stool collection tube (1) and use the stick to pick up sufficient sample quantity. Then, introduce the stick once into 4 different parts of the stool sample (2) making sure that, at each insertion, only the stick's screw is covered with the sample to take the appropriate amount of faecal sample and add it to the stool collection tube. Any additional sample that exceeds the stick's screw could cause wrong results. For liquid samples, add approx. 125µL in the stool collection tube using a micropipette.
2. Close the tube with the diluent and stool sample. Shake the tube in order to assure good sample dispersion (3).

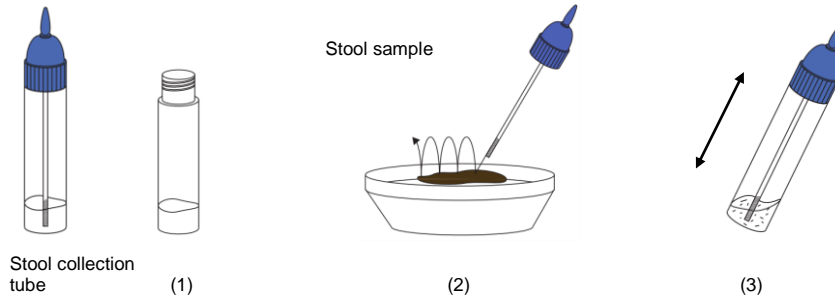**VERY IMPORTANT: Only fresh or fresh-frozen unpreserved and unfixed stool samples can be tested.**

### MATERIALS

#### MATERIALS PROVIDED

- CerTest *Crypto+Giardia+Entamoeba* combo card tests
- Instructions for use
- Stool collection tubes with diluent

#### MATERIALS REQUIRED BUT NO PROVIDED


- Specimen collection container
- Disposable gloves
- Timer

### TEST PROCEDURE

**Allow tests, stool samples and controls to reach room temperature (15-30°C) prior to testing. Do not open pouches until the performance of the assay.**

1. Proceed to shake the stool collection tube in order to assure good sample dispersion.
2. Remove CerTest *Crypto+Giardia+Entamoeba* combo card test from its sealed bag just before using it.
3. Take the stool collection tube, cut the end of the cap (4) and dispense 3 drops in the circular window marked with the letter A (5), 3 drops, using the same tube, in the circular window marked with the letter B (6), and 3 drops, using the same tube, in the circular window marked with the letter C (7). Avoid adding solid particles with the liquid.
4. **Read the results at 10 minutes.** Do not read the test result later than 10 minutes.

If the test does not run due to solid particles, stir the sample added in the sample window (S) with the stick. If it doesn't work, dispense a drop of extraction diluent until seeing the liquid running through the reaction zone.

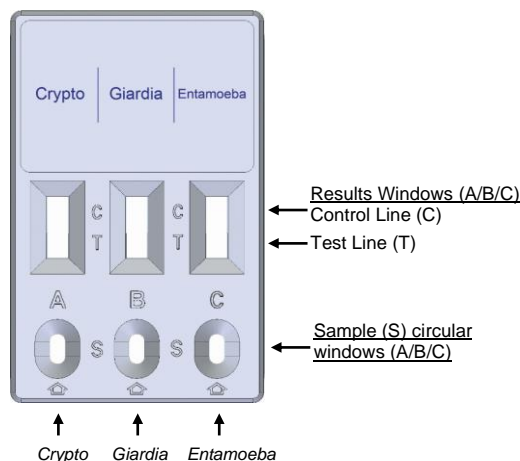
# CERTEST *Crypto*+*Giardia*+*Entamoeba*

One Step test to detect *Cryptosporidium*, *Giardia* and *Entamoeba* in combo card format /  
 Prueba combo de un solo paso para detección de *Cryptosporidium*, *Giardia* y *Entamoeba* en formato cassette


## *Crypto* Strip A-procedure

Add 3 drops in A-circular window

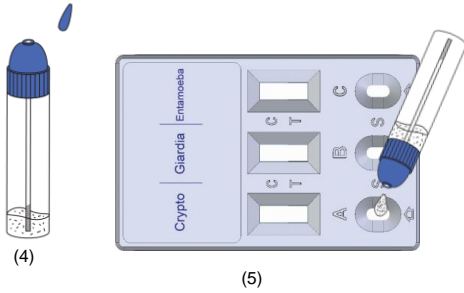
## *Giardia* Strip B-procedure

Add 3 drops in B-circular window

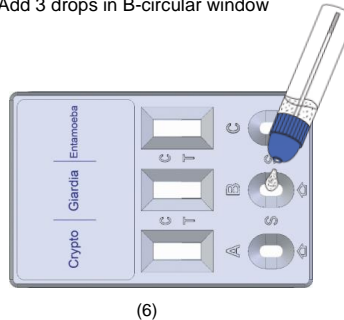
## *Entamoeba* Strip C-procedure


Add 3 drops in C-circular window

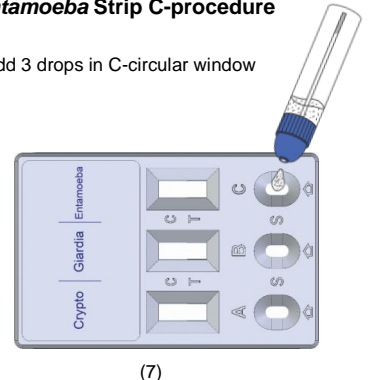
## INTERPRETATION OF THE RESULTS (please refer to the illustration below)

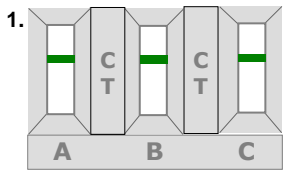
A: Green→Negative *Crypto*  
 B: Green→Negative *Giardia*  
 C: Green→Negative *Entamoeba*

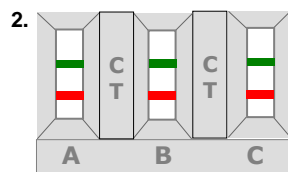
A: Green/Red→Positive *Crypto*  
 B: Green/Red→Positive *Giardia*  
 C: Green/Red→Positive *Entamoeba*

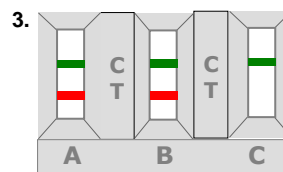
A: Green/Red→Positive *Crypto*  
 B: Green/Red→Positive *Giardia*  
 C: Green→Negative *Entamoeba*

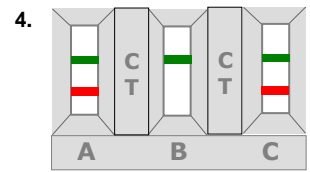
A: Green/Red→Positive *Crypto*  
 B: Green→Negative *Giardia*  
 C: Green/Red→Positive *Entamoeba*

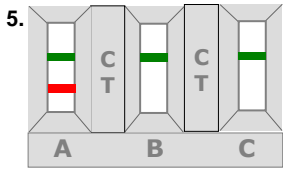
A: Green/Red→Positive *Crypto*  
 B: Green→Negative *Giardia*  
 C: Green→Negative *Entamoeba*

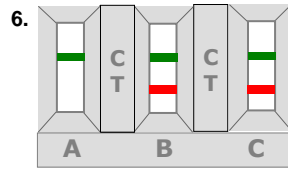A: Green→Negative *Crypto*  
 B: Green/Red→Positive *Giardia*  
 C: Green/Red→Positive *Entamoeba*
















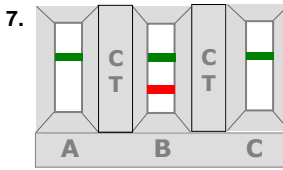
A: Green→Negative *Crypto*  
 B: Green/Red→Positive *Giardia*  
 C: Green→Negative *Entamoeba*

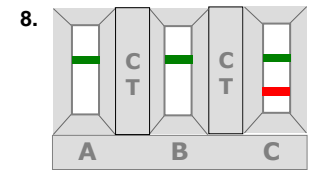
A: Green→Negative *Crypto*  
 B: Green→Negative *Giardia*  
 C: Green/Red→Positive *Entamoeba*

## INVALID RESULTS
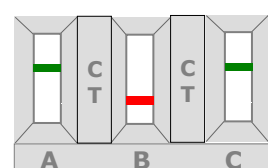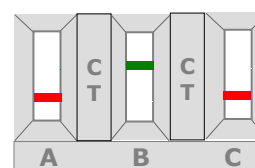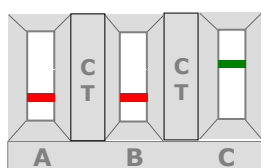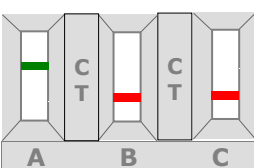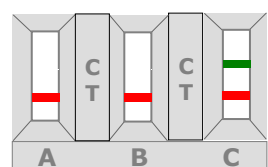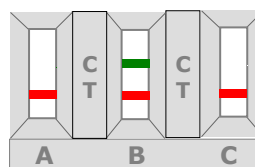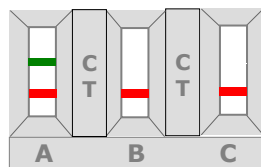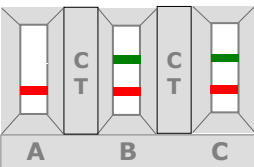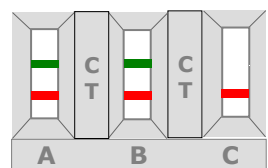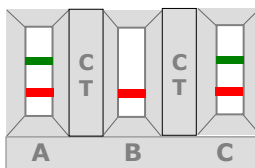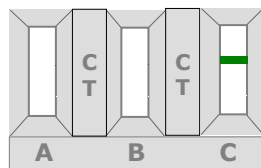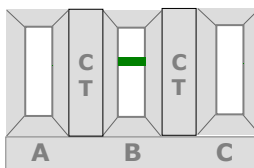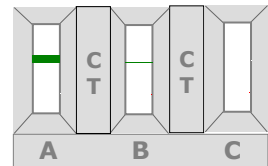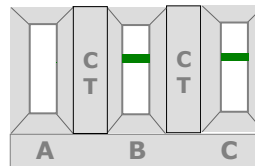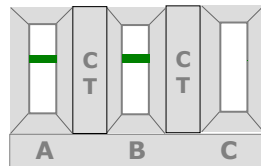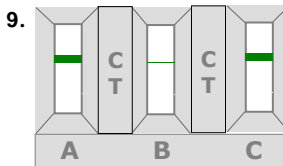
## CERTEST *Crypto*+*Giardia*+*Entamoeba*

One Step test to detect *Cryptosporidium*, *Giardia* and *Entamoeba* in combo card format/  
Prueba combo de un solo paso para detección de *Cryptosporidium*, *Giardia* y *Entamoeba* en formato cassette

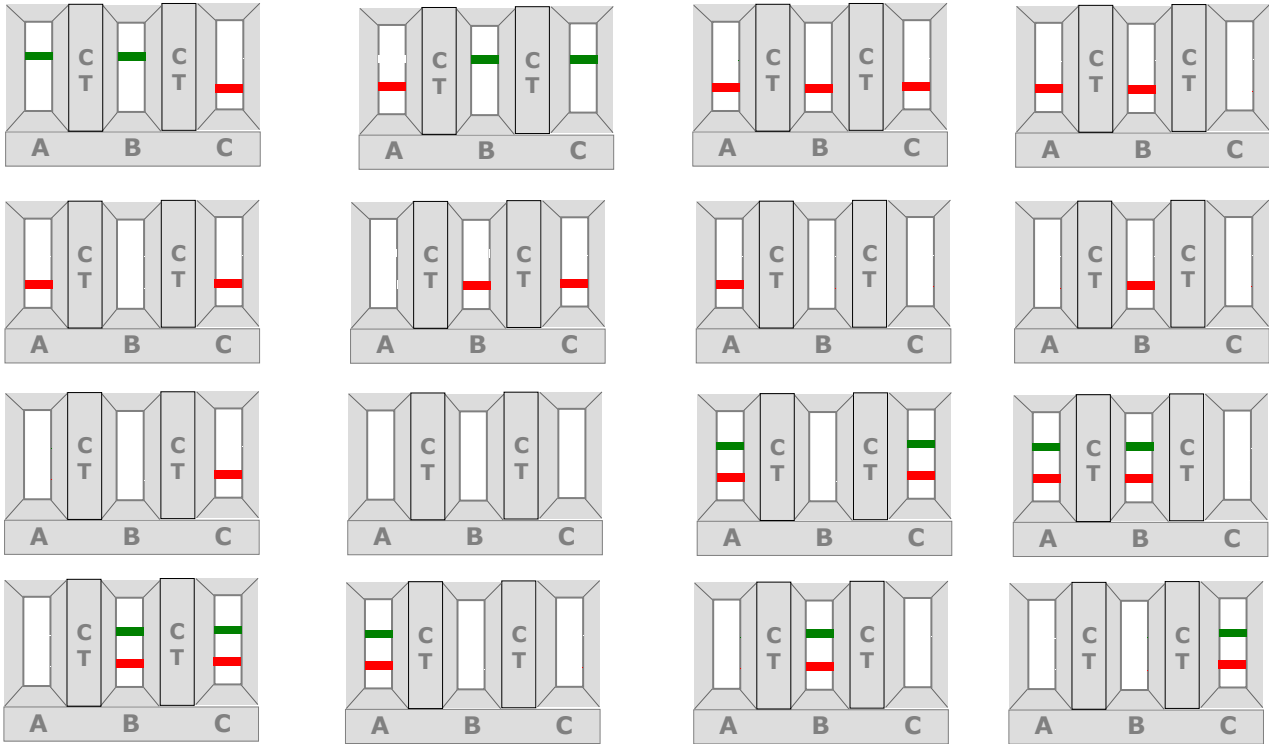	A ( <i>Crypto</i> )	B ( <i>Giardia</i> )	C ( <i>Entamoeba</i> )	Interpretation of the results
1.	- GREEN	- GREEN	- GREEN	There are no <i>Cryptosporidium</i> , <i>Giardia</i> or <i>Entamoeba</i> presence. No infection caused by <i>Cryptosporidium</i> , <i>Giardia</i> ( $\alpha$ 1- <i>giardin</i> and/or CWP1) or <i>Entamoeba</i> ( <i>histolytica</i> and/or <i>dispar</i> ).
2.	+ GREEN-RED	+ GREEN-RED	+ GREEN-RED	There are <i>Cryptosporidium</i> , <i>Giardia</i> and <i>Entamoeba</i> presence. Infection caused by <i>Cryptosporidium</i> , <i>Giardia</i> ( $\alpha$ 1- <i>giardin</i> and/or CWP1) and <i>Entamoeba</i> ( <i>histolytica</i> and/or <i>dispar</i> ).
3.	+ GREEN-RED	+ GREEN-RED	- GREEN	There is <i>Cryptosporidium</i> and <i>Giardia</i> presence. Infection caused by <i>Cryptosporidium</i> and <i>Giardia</i> ( $\alpha$ 1- <i>giardin</i> and/or CWP1).
4.	+ GREEN-RED	- GREEN	+ GREEN-RED	There is <i>Cryptosporidium</i> and <i>Entamoeba</i> presence. Infection caused by <i>Cryptosporidium</i> and <i>Entamoeba</i> ( <i>histolytica</i> and/or <i>dispar</i> ).
5.	+ GREEN-RED	- GREEN	- GREEN	There is <i>Cryptosporidium</i> presence. Infection caused by <i>Cryptosporidium</i> .
6.	- GREEN	+ GREEN-RED	+ GREEN-RED	There is <i>Giardia</i> and <i>Entamoeba</i> presence. Infection caused by <i>Giardia</i> ( $\alpha$ 1- <i>giardin</i> and/or CWP1) and <i>Entamoeba</i> ( <i>histolytica</i> and/or <i>dispar</i> ).
7.	- GREEN	+ GREEN-RED	- GREEN	There is <i>Giardia</i> presence. Infection caused by <i>Giardia</i> ( $\alpha$ 1- <i>giardin</i> and/or CWP1).
8.	- GREEN	- GREEN	+ GREEN-RED	There is <i>Entamoeba</i> presence. Infection caused by <i>Entamoeba</i> ( <i>histolytica</i> and/or <i>dispar</i> ).
9.	Any other result			Invalid result either A, B or C, we recommend repeating the assay using the same sample with another test.

**INVALID:** Total absence of any control coloured line (GREEN) regardless the appearance or not of the test lines (RED). Insufficient specimen volume, incorrect procedural techniques or deterioration of the reagents are mostly the main reasons for control lines failure. Review the procedure and repeat the assay with a new test. If the symptoms or situation still persist, discontinue using the test kit and contact your local distributor.

## CERTEST *Crypto+Giardia+Entamoeba*

One Step test to detect *Cryptosporidium*, *Giardia* and *Entamoeba* in combo card format /  
Prueba combo de un solo paso para detección de *Cryptosporidium*, *Giardia* y *Entamoeba* en formato cassette

### NOTES ON THE INTERPRETATION OF RESULTS

The intensity of the red coloured bands in the test lines (T) in the results windows will vary depending on the concentration of antigens present in the specimen. However, neither the quantitative value nor the rate of increase in antigens can be determined by this qualitative test.

### QUALITY CONTROL

Internal procedural controls are included in the test. The green lines appearing in the control lines (C) in the results windows are internal controls, which confirm sufficient specimen volume and correct procedural technique.

### LIMITATIONS

1. The test must be carried out within 2 hours after opening the sealed bag.
2. Only fresh or fresh-frozen unpreserved and unfixed stool samples can be tested.
3. An excess of sample could cause wrong results (brown bands appear). Dilute the sample with the diluent and repeat the test.
4. The intensity of test line may vary from very strong at high antigens concentration to faint when the antigens concentration is close to the detection limit value of the test.
5. After one week of infection, the number of parasites in faeces is decreasing, making the sample less reactive. Stool samples should be collected within one week of the onset of symptoms.
6. CerTest *Crypto+Giardia+Entamoeba* should be used only with samples from human faeces. The use of other samples has not been established. The quality of the test depends on the quality of the sample; proper faecal specimens must be obtained.
7. Positive results determine the presence of *Cryptosporidium*, *Giardia* ( $\alpha$ 1-*giardin* and/or CWP1), *Entamoeba* (*histolytica* and/or *E. dispar*) in faecal samples, nevertheless, a positive result should be followed up with additional laboratory techniques (biochemical methods or microscopy) to confirm the results. A confirmed infection should only be made by a physician after all clinical and laboratory findings have been evaluated and must be based in the correlation of the results with further clinical observations.
8. A negative result is not meaningful because of it is possible the antigens concentration in the stool samples is lower than the detection limit value. If the symptoms or situation still persist, *Cryptosporidium*, *Giardia*, *Entamoeba* determination should be carried out with another technique (for example microscopy).
9. Bloody stool samples can contain components that may cause non-specific reactions in the test. Every bloody stool sample whose result is positive should be followed up with other techniques of diagnosis to confirm the result.
10. Mucous stool samples may cause non-specific reactions in the test- Every mucous stool sample whose result is positive should be followed up with other techniques of diagnosis to confirm the result.

### EXPECTED VALUES

Parasitic diseases are incriminated in causing more than 33% of global deaths of which intestinal parasitic infections are believed to take the major share. Lack of safe drinking water and environmental sanitation are largely responsible for more than 800 million expected cases of diarrheal diseases and 4.5 million associated deaths in many developing countries every year (3).

Morbidity and mortality due to diarrheal diseases in developing countries remain to be the main public health problems that need due attention. Although there could be many other causes of diarrhoea, the enteric protozoa *Cryptosporidium parvum* and *Giardia lamblia* have been recognized as important causes of both out break-related and sporadic diarrhoea among human beings. Both immunocompetent and immunocompromised individuals could be the victims of diarrheal diseases caused by these parasites.

*Entamoeba histolytica* is the third main cause of morbidity and mortality due to parasitic disease in humans after malaria and schistosomiasis. It is estimated that *E. histolytica* and *E. dispar* parasitize approximately 10% of the world population, of which 90% are asymptomatic infections. About between 50000 and 100000 deaths every year are due to this disease.

In developed countries, amoebiasis tends to be more common in older and immunodepressed individuals. However, in tropical regions, the epidemiology of amoebiasis is completely different and is more common among the general population and particularly among patients attending health care centres with diarrhoea.

Out of 10% of the world's population infected by *E. histolytica*, only 1% develops invasive form of the disease. Only approximately 10% (pathogenic infections) need to be treated. In developed countries, *E. histolytica* infection is not endemic and occurs mostly after traveling to areas of endemicity.

### PERFORMANCE CHARACTERISTICS

#### Clinical sensitivity and specificity

An evaluation, with faecal samples, was performed comparing the results obtained by an immunochromatographic test (CerTest *Crypto+Giardia+Entamoeba*, CerTest) vs the Evaluation Criteria. Evaluation criteria: two rapid test test were evaluated (CerTest *Crypto+Giardia+Entamoeba*; CerTest (Crypto strip) and other rapid test from the competitor); the discrepant results were confirmed by qPCR technique (Viasure *Cryptosporidium, Giardia & E. histolytica* Real Time Detection Kits (CerTest). The results were as follows:

		Evaluation Criteria (Crypto)		
		+	-	Total
IC test: CerTest <i>Crypto+Giardia+Entamoeba</i> (Crypto)	+	33	0	33
	-	2	91	93
	Total	35	91	126

CerTest <i>Crypto+Giardia+Entamoeba</i> (Crypto) vs Evaluation Criteria		
	Mean Value	95% CI (confidence interval)
<b>Sensitivity</b>	94.3%	80.8 – 99.3%
<b>Specificity</b>	100.0%	96.0 – 100.0%
<b>PPV</b>	100.0%	89.4 – 100.0%
<b>NPV</b>	97.8%	92.4 – 99.7%

An evaluation, with faecal samples, was performed comparing the results obtained by an immunochromatographic test (CerTest *Crypto+Giardia+Entamoeba*, CerTest) vs the Evaluation Criteria. Evaluation criteria: three rapid test were evaluated (CerTest *Crypto+Giardia+Entamoeba*, Certest (Giardia strip) and other two rapid test from the competitors); the discrepant results were confirmed by qPCR technique (Viasure *Cryptosporidium, Giardia & E. histolytica* Real Time Detection Kits (CerTest). The results were as follows:

		Evaluation Criteria (Giardia)		
		+	-	Total
IC test: CerTest <i>Crypto+Giardia+Entamoeba</i> (Giardia)	+	44	2	46
	-	1	79	80
	Total	45	81	126

CerTest <i>Crypto+Giardia+Entamoeba</i> (Giardia) vs Evaluation Criteria		
	Mean Value	95% confidence interval
<b>Sensitivity</b>	97.8%	88.2–99.9%
<b>Specificity</b>	97.5%	91.4–99.7%
<b>PPV</b>	95.7%	85.2–99.5%
<b>NPV</b>	98.8%	93.2–100.0%

## CERTEST *Crypto+Giardia+Entamoeba*

One Step test to detect *Cryptosporidium*, *Giardia* and *Entamoeba* in combo card format/  
Prueba combo de un solo paso para detección de *Cryptosporidium*, *Giardia* y *Entamoeba* en formato cassette

An evaluation, with faecal samples, was performed comparing the results obtained by an immunochromatographic test (CerTest *Crypto+Giardia+Entamoeba*, CerTest) vs the Evaluation Criteria. Evaluation criteria: two rapid test were evaluated (CerTest *Crypto+Giardia+Entamoeba*, CerTest (*Entamoeba* strip) and other rapid test from the competitors); the discrepant results were confirmed by qPCR technique (Viasure *Cryptosporidium*, *Giardia* & *E. histolytica* Real Time Detection Kits and Viasure *Entamoeba dispar* Real Time Detection kit (CerTest). The results were as follows:

	Evaluation Criteria ( <i>Entamoeba</i> )			
		+	-	Total
IC test: CerTest <i>Crypto+Giardia+Entamoeba</i> ( <i>Entamoeba</i> )	+	10	0	10
	-	3	65	68
	Total	13	65	78

CerTest <i>Crypto+Giardia+Entamoeba</i> ( <i>Entamoeba</i> ) vs Evaluation Criteria		
	Mean Value	95% confidence interval
<b>Sensitivity</b>	76.9%	46.2–95.0%
<b>Specificity</b>	100.0%	94.5–100.0%
<b>PPV</b>	100.0%	69.2–100.0%
<b>NPV</b>	95.6%	87.6–99.1%

The results showed a high sensitivity and specificity to detect *Cryptosporidium*, *Giardia* ( $\alpha 1$ -*giardin* and/or CWP1) and *Entamoeba* (*histolytica* and/or *dispar*) using CerTest *Crypto+Giardia+Entamoeba*.

### Cross reactivity

An evaluation was performed to determine the cross reactivity of CerTest *Crypto+Giardia+Entamoeba*; no cross reactivity against gastrointestinal pathogens, other organism, substances and/or faecal markers occasionally present in faeces:

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- MARSHALL, M.M., et al., "Waterborne Protozoan Pathogens", *Clinical Microbiology Review*, Jan. 1997, pp 67-85.
- DYLAN R. PILLAI and KEVIN C, KAIN, "Immunochromatographic Strip-Based Detection of *Entamoeba histolytica*-*E. dispar* and *Giardia lamblia* Coproantigen". *Journal of Clinical Microbiology*, Sept. 1999, Vol. 37, No 9, p. 3017-3019.
- MARTIN-AMPUDIA M., Under-notification of cryptosporidiosis by routine clinical and laboratory practices among non-hospitalised children with acute diarrhoea in Southern Spain", *Journal of Infection*, 2012, p. 40:113-119.
- LYNNE S. GARCIA et al., "Commercial Assay for Detection of *Giardia lamblia* and *Cryptosporidium parvum* Antigens in Human Fecal Specimens by Rapid Solid-Phase Qualitative Immunochromatography", *Journal of Clinical Microbiology*, Jan. 2003, Vol. 41, No. 1, p. 209-212.
- BLESSMANN J., et al. "Longitudinal Study of Intestinal *Entamoeba histolytica* Infections in Asymptomatic Adult Carriers". *Journal of clinical microbiology*, Oct. 2003, p. 4745-4750 Vol. 41, No. 10.
- GONIN P., et al. "Detection and Differentiation of *Entamoeba histolytica* and *Entamoeba dispar* Isolates in Clinical Samples by PCR and Enzyme-Linked Immunosorbent Assay". *Journal of clinical microbiology*, Jan. 2003, p. 237-241 Vol. 41, No. 1.
- ADAM RODNEY. Biology of *Giardia lamblia*. *Clin Microbiol Rev*. 2001;14(3):447-475. doi:10.1128/CMR.14.3.447-475.2001).
- GAËTAN FAUBERT. Immune Response to *Giardia duodenalis*. *Clin Microbiol Rev*. 2000; 13(1): 35-54. 0893-8512/00/\$04.00+.

## SYMBOLS FOR IVD COMPONENTS AND REAGENTS/SÍMBOLOS PARA REACTIVOS Y PRODUCTOS PARA DIAGNÓSTICO IN VITRO

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